

User Manual

Python Recording Aid

Contents

1	Introduction	2
2	Installation	2
2.1	Installing anaconda	2
2.2	For Windows users:	2
2.3	For Linux users:	2
2.4	yml file installation:	3
3	Main Menu	3
3.1	Loading Image Folder	3
3.2	Load Series	3
3.3	Enabling Extra Options	3
3.4	Reslice	4
3.5	Add Electrode(Chamber)	4
3.6	Flip Menu	4
4	Coregistration Menu	4
4.1	Flip and Swap options	4
4.2	Send to Main Form	6
5	Some usable notes	6

1 Introduction

Imaging guided access to the brain has become a routine procedure for various research and clinical applications, including drug administration, neurophysiological recording, and sampling tissue. Python Recording Aid (PyRecAid) is an open-source tool designed to automate accessing the area of interest in brain.

The software supports both Dicom and Nifti image format as input and provides processing procedures including reslicing, flipping, adding chamber, and coregistration in an automated manner.

2 Installation

Simply download or clone the directory to a desired location on your local system (extract the .zip file in case of downloading). Then redirect to this folder on your system to find the following files and folders:

src → here you will see the program scripts including the main method [...]. files → here you will see the manual.pdf containing user guide.

Prior to the installation you should have a distribution of Python3 interpreter installed on your system (we recommend the anaconda distribution). Instructions on how to run the program methods on Windows and Linux machines are described below. In the case you do not want to install the packages manually, after installing anaconda, jump to the yml installation part.

2.1 Installing anaconda

Download the anaconda from [here](#) for your operating system. Then install then application and follow bellow procedure to install software requirements.

2.2 For Windows users:

After installation of Python3, open command prompt or the anaconda prompt (depending on your Python3 distro and path settings) and type the following commands at the prompt (case sensitive):

```
$ pip install scipy dicom PyQt5 matplotlib OpenCV-Python scikit-image  
transforms3d dill nibabel torch torchvision torchaudio
```

If you have python2 in your environment then you should type pip3 instead of pip in the above command since in that case pip defaults to Python2. It takes some time so be patient and wait for the transactions to finish and then enter:

```
$ conda install -c clinicalgraphics gdc
```

In case you do not want to install anaconda, you can download the GDCM library at [here](#).

2.3 For Linux users:

Linux machines contain a system Python3. It is highly recommended that you don't install the program using the system Python3. You are better to install a virtual environment by venv or conda (in case you have anaconda or miniconda) and then follow the instructions below. Open a terminal with your desired python3 environment activated and type in the following line:

```
$ pip install scipy==1.2.0 dicom PyQt5 matplotlib OpenCV-Python scikit-image  
transforms3d dill nibabel torch==1.10.2+cpu torchvision==0.11.3+cpu  
torchaudio==0.10.2+cpu -f  
https://download.pytorch.org/whl/cpu/torch_stable.html
```

If you have python2 in your environment then you should type pip3 instead of pip in the above command since in that case pip defaults to Python2. It takes some time so be patient and wait for the transactions to finish and then enter:

```
$ conda install -c clinicalgraphics gdc
```

In case you do not want to install anaconda, you can download the GDCM library at [here](#).

2.4 yml file installation:

Installing from yml file is the easiest way to execute the software. After installing the anaconda open a terminal (Linux/Mac) or anaconda command prompt (Windows) and create a new environment by running:

```
$ conda env create -f environment.yml
```

By running this command a new environment will be created under the name "pyRecAid". By executing the bellow command you can see the existing environments in your device.

```
$ conda env list
```

If the installation process runs successfully, you should be able to see the "pyRecAid" environment under the exported list. After installing successfully, you have to activate the environment by running:

```
$ conda activate pyRecAid
```

After activating the environment, go to the cloned code directory on your system and run:

```
$ python3 Main.py
```

Thereafter, the software starts and ready to be used.

3 Main Menu

This menu contain tools to show the 3d images, reslice it, add chamber, and flip it if required. In the following sections we express each part by details:

3.1 Loading Image Folder

Using **Load Folder** button, the series directory is selected. If there exist image series inside the selected directory, the list of existing series will be showed as shown in figure 1.

3.2 Load Series

After loading all images inside directory, you need to select the desired series and then press **Load Series** button. Then the loaded series is showed in main window as figure 2. As it is observable, images contain the the slice distances in milimeters. Three slices including Horizontal, Sagittal, and Coronal are showed. By moving the mouse scroll on each slice, the whole brain could be covered.

3.3 Enabling Extra Options

To enable reslice, add chamber, and flip the tools menu should be opened from the upper bar of frame. Then, by pressing each option the related window will be shown on the side bar of main window (figure 3)

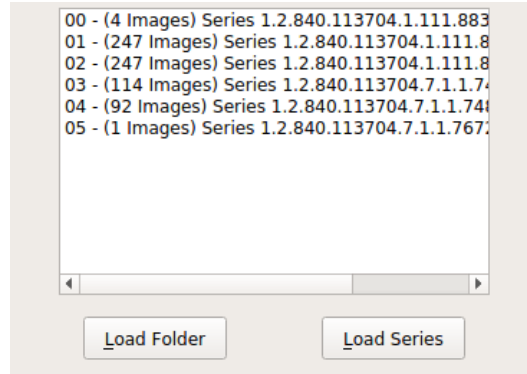


Figure 1: Series List

3.4 Reslice

This tool is used to reslice the brain by specifying AP and ML angles. Both angles ranges are from -180 to 180 degrees. To reslice the brain, both angles should be determined first, then by pressing **Reslice** button a new window with resliced series is shown as figure 4.

3.5 Add Electrode(Chamber)

Adding an electrode to the brain is one of the options that is provided in both resliced and main menu. As shown in figure 5, an electrode is added in the reslice window by setting required parameters. These parameters could be defined in the **Electrode Size (Chamber Size)** window (figure 6).

3.6 Flip Menu

This menu is very handy when the slices are not shown in right Directions. If the directions of head in series are not as desired, then by flipping them they could be put in right directions. As a case in point, Sagittal view is flipped as shown in figure 7.

4 Coregistration Menu

This menu is accessible from both main window by pressing **Coregistration** button, or opening from tools bar. The coregistration window is shown in figure 8. As it is observable, there are three rows for three series. The first and second series are the two series we intend to register. Consequently, the third row is for result series, which is the combination both first and second series. Loading the series is quite similar to the main window. After loading both series, four match points are required to register both series. As the utilized algorithm is based on 3D affine transformation, at least one of the four match points should be on a different slice (not all of four points on a same plane). Point selection window is designed for this purpose as shown in figure 9. Points could be selected by setting the coordinates on both images by mouse scroll and the pressing **Select Point** button. By selecting each point, a check mark will appear in each box related to the point. If the user wants to clear points, it could be done by pressing **Clear Points** button. If a single point is to be cleared, then the user should clear the check mark of corresponding point.

4.1 Flip and Swap options

If the series are not in the same directions, then the user could change it to the right order by using **Flip** or **Swap** buttons. Before using this tools, the related series should be chosen and then the transformation could be applied

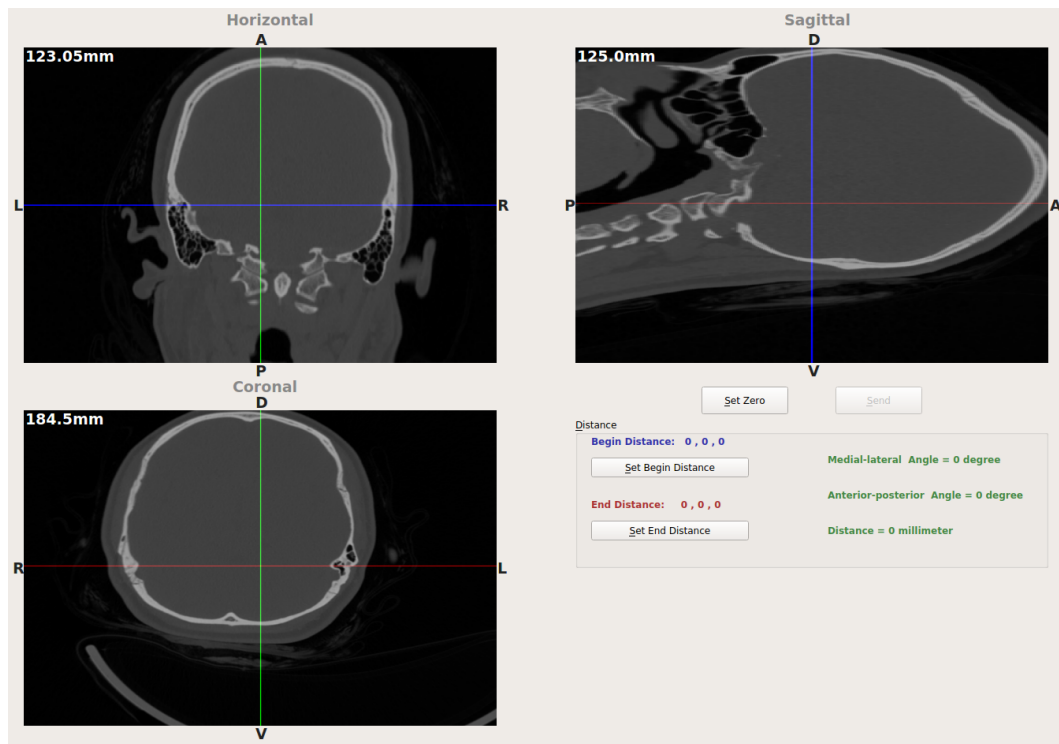


Figure 2: Loaded Series

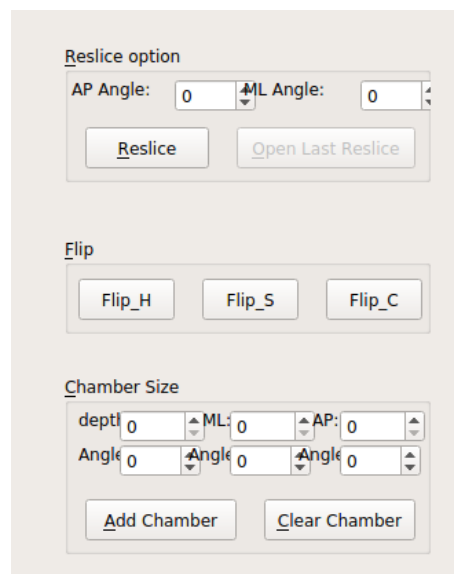


Figure 3: Extra Options

(figure 10). These option could be enabled or disabled from tools bar. Two series should be in a same direction in order to reach a suitable coregistration.

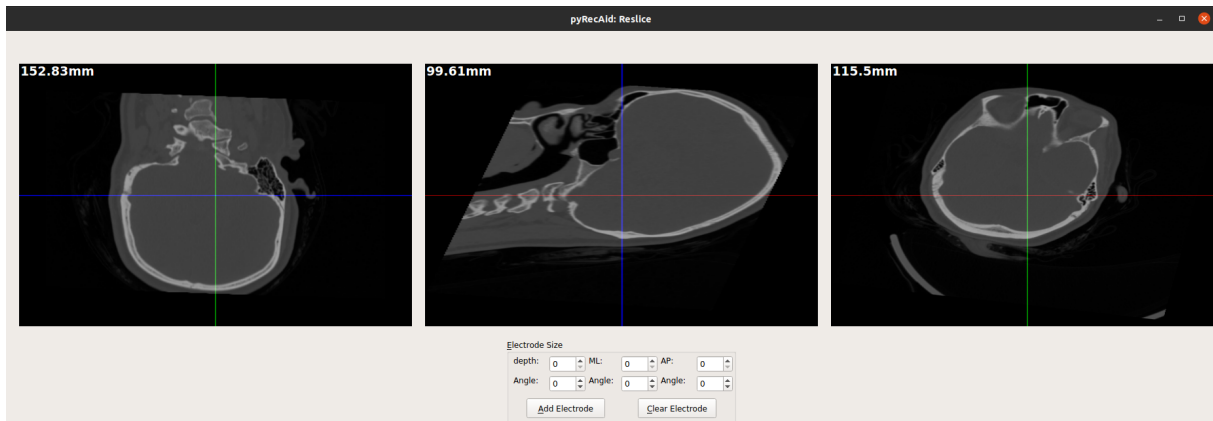


Figure 4: Reslice Window

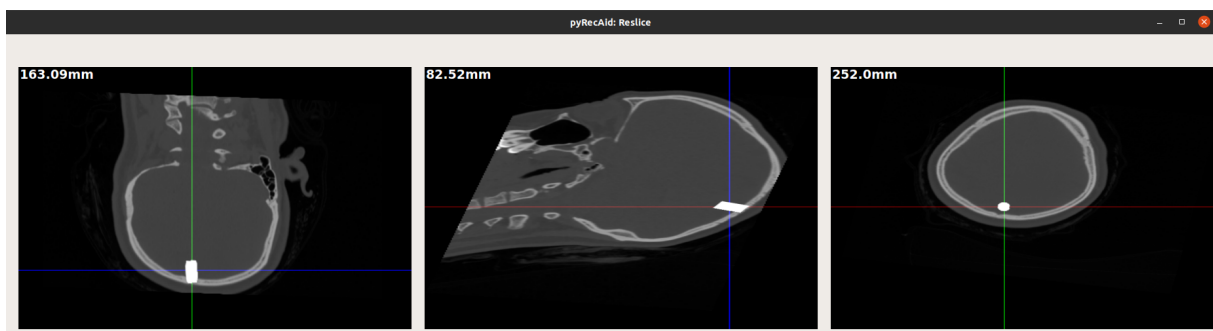


Figure 5: Adding Electrode

Electrode Size

depth:	10	ML:	10	AP:	40
Angle:	5	Angle:	25	Angle:	5


Figure 6: Electrode Size

4.2 Send to Main Form


After applying coregistration on series, the result could be sent to the main window by pressing **Send To Main Form** button. In the main window, the series are available for reslicing or adding chamber.

5 Some usable notes

- Each slice could be maximized by double clicking on it.

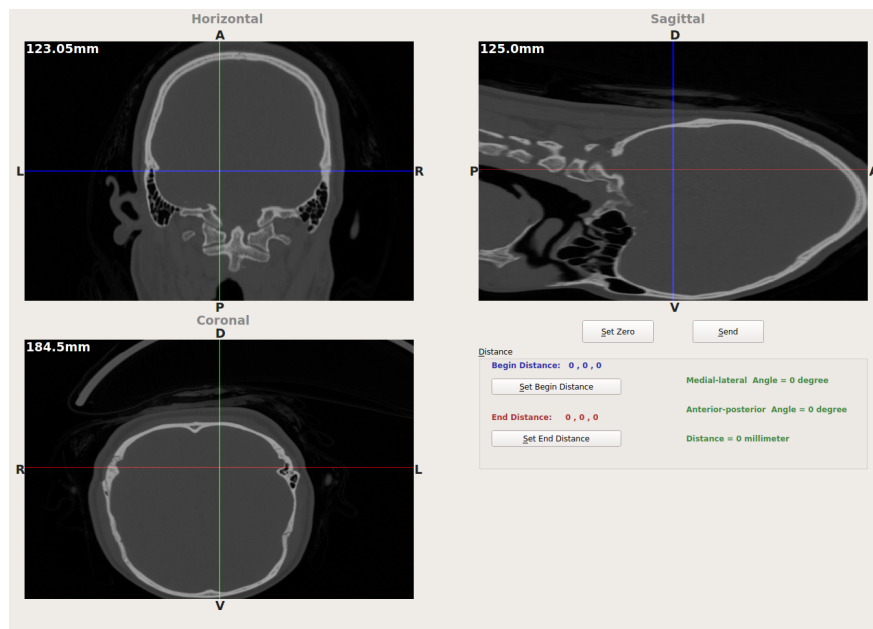

(a)

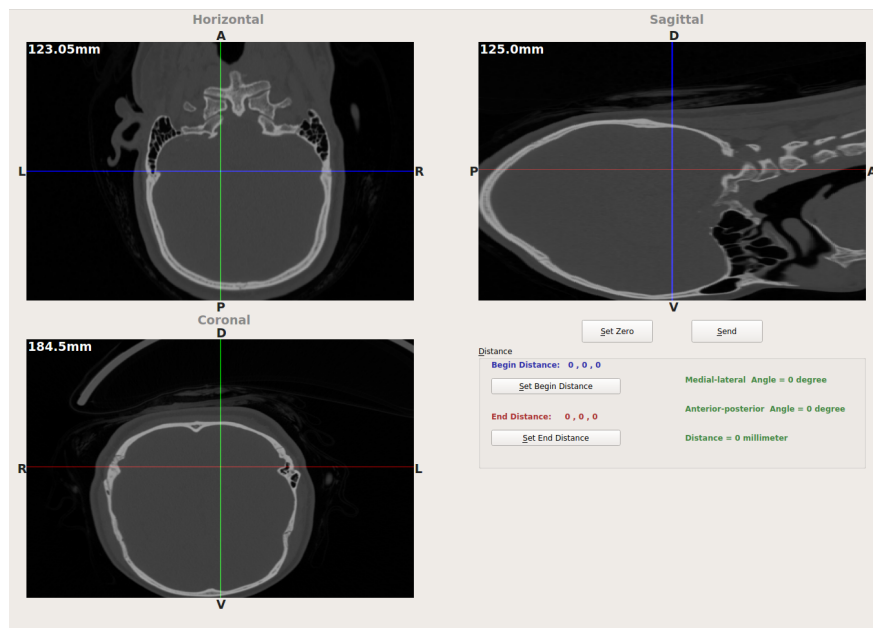

(b)

Figure 7: Flipped images

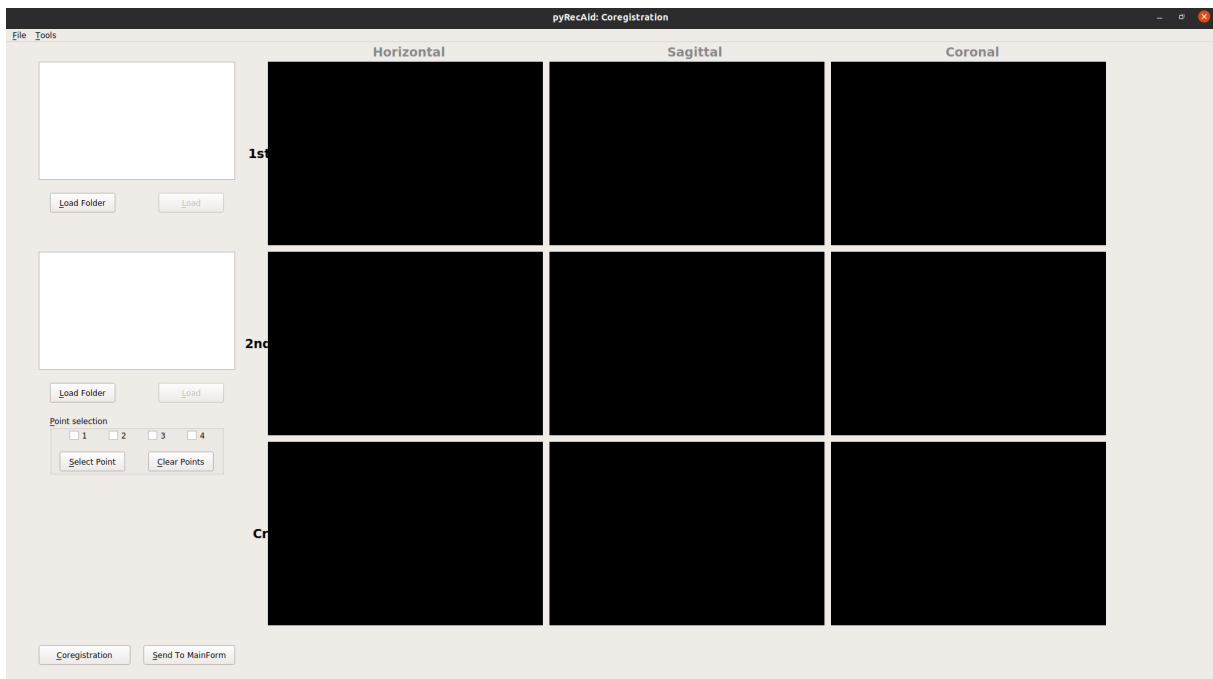


Figure 8: Coregistration Menu

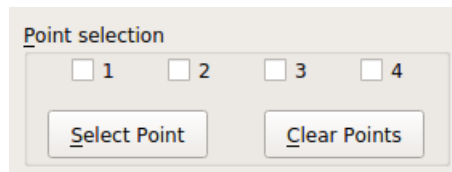


Figure 9: Point Selection

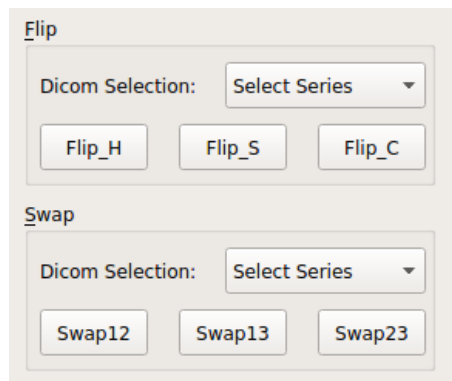


Figure 10: Point Selection